




Paraguard Excel Stretcher

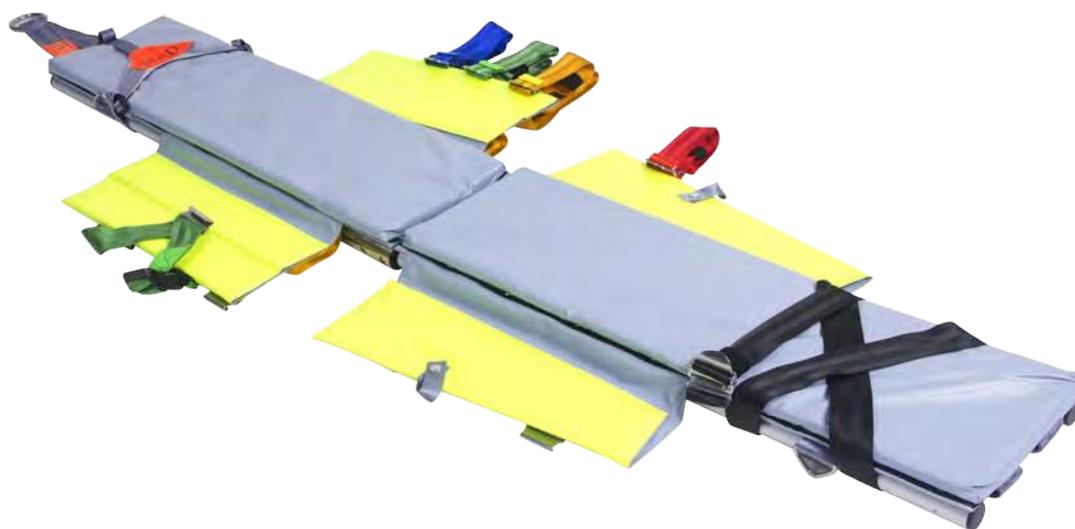

The Paraguard Stretcher is a compact and versatile rescue stretcher, designed for a variety of emergencies. It is ideal for search and rescue operations and can be used to rescue casualties from building collapse, industrial incidents or confined spaces.

The Paraguard Excel Stretcher has been developed from the popular Paraguard MKIII to provide a superior rescue stretcher. Constructed from flame resistant materials, high quality aluminum and stainless steel, the Paraguard Stretcher is impervious to water, grease, oil and other petroleum products and is resistant to corrosion and rot.

It is particularly useful for high building, mountain and helicopter rescues. A 4-Point Lifting harness is available with sturdy yet lightweight karabiners for vertical lifting. In use since WWII Paraguard medical stretchers are approved for use by many authorities including NATO 6530-99-446-9403.

The Paraguard Stretcher is quick and simple to use with various straps to secure the chest, arms, thighs and lower legs. The head is held in place by a 3-point, non-slip forehead strap. A figure of eight strap supports the feet. To enhance infection control the flaps are wipe clean and the patient pads may be detached from the base for washing.

Also available for the Paraguard Stretcher are a set of slot in carrying handles so that two medics can carry the victim. These can also be used in conjunction with the Paraguard Shoulder Harness which will displace the victim's weight onto the rescuers shoulders, hips and legs.

When moving in awkward spaces the unique construction allows the centre of the Paraguard stretcher to unlock and bend around corners. For fast deployment new speed links also enable the stretcher to quickly unfold when removed from storage. For storage and ease of transport this medical stretcher folds in half to fit into a back-pack (aka valise) supplied as standard.

Golden Season Pte Ltd



Paraguard Excel Stretcher - Specification



Specifications	
Length:	182cm (92cm folded)
Width:	25cm (27cm folded)
Depth:	15cm (Folded only)
Net Weight:	11.5kg
Max. Load:	136kg

*Specifications are subject to change without prior notice.